

Rapid Malaria Pf (HRP 2)/ PAN (pLDH) Antigen Test - Device

For *In-Vitro* Diagnostic Use Only

Store at 4°C to 40°C

OVERVIEW

Malaria is a serious parasitic disease characterized by fever, chills, and anaemia and is caused by a parasite that is transmitted by the bite of infected Anopheles mosquitoes. There are four kinds of malaria parasites that can infect humans: *Plasmodium falciparum*, *P. vivax*, *P. ovale*, and *P. malariae*. In humans, the parasites (called sporozoites) migrate to the liver where they mature and release another form, the merozoites. At present Malaria is diagnosed microscopically using thick and thin blood films. These require expert knowledge to correctly identify the species, not always available 24 hrs a day. This is where a reliable support test becomes invaluable.

INTENDED USE

This is a rapid, in vitro, qualitative lateral flow immunoassay for the detection of *P. falciparum* specific Histidine rich protein-2 (Pf, HRP-2) and PAN malaria specific pLDH from human whole blood samples. The test may also be used for the differentiation of *P. falciparum* and other malarial species and for the follow up of anti-malarial therapy in whole blood samples.

PRINCIPLE

The Rapid Malaria Pf (HRP2)/ PAN (pLDH) Antigen test contains a membrane strip, which is pre-coated with two test lines and one control line. One monoclonal antibody (Pf - test line 1) is specific to Pf Histidine Rich Protein 2 (HRP2) of the *Plasmodium falciparum* species and the other line (PAN - test line 2) consists of a monoclonal antibody specific to PAN plasmodium lactate Dehydrogenase (pLDH). The control line (C) consists of Goat anti-Rabbit IgG. The conjugate pad is dispensed with HAMA blocking reagent and colloidal gold conjugated to *P. falciparum* (HRP2) specific, PAN (pLDH) specific, antibodies and rabbit IgG. The test is designed for the differential diagnosis between *Plasmodium falciparum* and the other Pan specific species (*P. vivax*, *malariae* and *ovale*).

After addition of the blood sample and the assay buffer to the respective wells on the test containing a test strip, the whole blood gets lysed and if the sample contains detectable levels of the Pf (HRP2) antigen and/or PAN (pLDH) antigen it reacts with the respective gold conjugated with malaria Pf (HRP2) specific antibodies and/or PAN (pLDH) specific antibodies to form a complex. The unbound colloidal gold particles along with complex move on to the nitrocellulose membrane. This complex moves further and reacts with the respective malaria Pf (HRP2) specific antibodies/ PAN (pLDH) specific antibodies test lines on the nitrocellulose membrane area to form a colored band (Test band/s). The unbound complex, unbound gold and the rabbit IgG conjugated colloidal gold particles move further to the goat-anti rabbit IgG coated control area to form a colored band (C- Control line). The appearance of test lines and control line in respective area indicates the positive result. Appearance of only control line indicates a negative result. The control line acts as a procedural control. Control line should always appear if the test is performed as per the procedure and reagents are working properly.

MATERIAL PROVIDED

1. Test: Nitrocellulose Membrane assembly pre-dispensed with monoclonal anti-Pf (HRP-2) antibody, monoclonal PAN (pLDH) antibody, Goat anti rabbit IgG, Conjugate strip containing HAMA blocking reagent and colloidal gold conjugated monoclonal anti-Pf (HRP-2) antibody, PAN (pLDH) antibody, and rabbit IgG at the respective regions.
2. Desiccant pouch
3. Disposable 5µl sample applicator
4. Package Insert
5. Assay Buffer

OPTIONAL MATERIAL REQUIRED

1. Calibrated micropipette capable of delivering 5 µl sample accurately.
2. Stop watch.
3. Disposable gloves

PRECAUTIONS/KIT STORAGE AND STABILITY

1. Please read all the information in this package insert before performing the test. Pay particular attention to the position of the Control and Test lines.
2. Do not use after the expiration date printed on the foil pouch.
3. Store in the sealed pouch in a dry place in between temperature 4°C to 40°C. Do not freeze.
4. Do not use if pouch is torn or damaged.
5. Do not open the foil pouch until you are ready to start the test.
6. Keep out of the reach of children.

WARNINGS

1. Do not reuse the test.
2. Follow the instruction to get accurate results.
3. Use appropriate personal protective equipment.
4. Dispose of hygienically as per local regulatory requirements.
5. Do not touch the membrane.
6. Treat blood samples and used tests as potentially infectious. Avoid contact with skin.
7. For in vitro diagnostic use. Not to be taken internally.
8. Do not eat the desiccant in the package.
9. Do not mix the specimen sample or interchange the different specimen.

SPECIMEN COLLECTION

Fresh anti coagulated whole blood should be used as a test sample. EDTA or Heparin or Oxalate or Tri-sodium Citrate can be used as suitable anticoagulants. The specimen should be collected in a clean glass or plastic container. If immediate testing is not possible then store the specimen at 2°C to 8°C for up to three days before testing. Clotted or contaminated blood samples should not be used for performing the test. Fresh blood from finger prick/ puncture may also be used as a test specimen.

TEST PROCEDURE

1. Bring the kit components to room temperature before testing.
2. Open the pouch and retrieve the test and desiccant pouch. Check the color of the desiccant. It should be blue, if it has turned colorless or pink, discard the test and use another test. Once opened, the test must be used immediately.
3. Label the test with patient's identity.
4. Tighten the vial cap of the assay buffer provided with the kit in the clockwise direction to pierce the dropper bottle nozzle.
5. Evenly mix the anti-coagulated blood sample by gentle swirling. Dip the sample loop into the sample. Ensuring that a loop full of blood is retrieved, blot the blood so collected in the sample port 'S'. (This delivers approximately 5µl of the whole blood specimen).

OR

In case finger prick blood is being used, touch the sample loop to the blood on the finger prick. Ensuring that a loop full of blood is retrieved, immediately blot the specimen in the sample port 'S'. (Care should be taken that the blood sample has not clotted and the transfer to the sample port is immediate).

OR

Alternatively, 5µl of the anti-coagulated or finger prick specimen may be delivered in the sample port 'S' using a micro pipette.


NOTE: Ensure that the blood from the sample loop has been completely taken up at the sample port 'S'.

6. Immediately dispense three drops of assay buffer (Approx. 90 µl) into buffer port 'B', by holding the plastic dropper bottle vertically.

7. Read the results at the end of 20 minutes. Do not read the result after 30 minutes.

INTERPRETATION OF RESULTS

NEGATIVE for Malaria : If colored band appears at the control region 'C'.

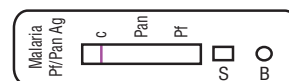


POSITIVE for *P. falciparum* and/or mixed infection: In addition to the control band, two Colored bands appear at regions 'Pf' and 'Pan' in the test window.



POSITIVE for Other species (non falciparum): In addition to the control band, one Colored band appears only at region 'Pan' in the test window.




POSITIVE for *P. falciparum*: In addition to the control band, one Colored band appears only at region 'Pf' in the test window.

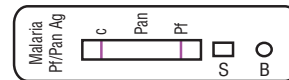

INVALID : The test should be considered invalid if,


A) No line appears at 'C', 'Pf' and 'Pan' regions.

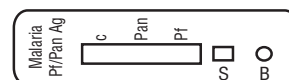


B) No line appears at 'C' region and line appear at 'Pf' and 'Pan' region.



C) No line appears at 'C' and 'Pf' region and line appear at 'Pan' region.

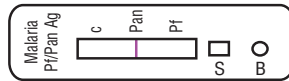

D) No line appears at 'C' and 'Pan' region and line appear at 'Pf' region.

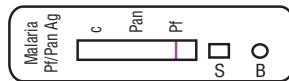

Performance Characteristics

Internal Evaluation :

In an in-house study, total 200 samples were evaluated for sensitivity and specificity. We found the relative sensitivity was 100 % (i. e. 100/100) and the relative specificity was 100 % (i. e. 100/100). The results are summarized in the following table :

Sample	Total Number of Samples Tested	Rapid Malaria Pf (HRP 2)/PAN (pLDH) Antigen Test		Sensitivity %	Specificity %
		Positive	Negative		
Malaria Pf Positive Whole Blood Samples	52	52	0	100	-
Malaria Pan Positive Whole Blood Samples	68	68	0	100	-
Malaria Negative Whole Blood Samples	100	0	100	-	100

Cross reactivity was studied using RF positive samples and no cross reactivity was observed.

External Evaluation:

In an external study, total 200 samples were evaluated for sensitivity and specificity. Relative sensitivity was 100% (i. e. 50/50) and the relative specificity was 100 % (i. e. 150/150). Positive Predictive Value (PPV) and Negative Predictive Value (NPV) for the test was 100%. The results are summarized in the following table:

Sample	Total Number of Samples Tested	Rapid Malaria Pf (HRP 2)/PAN (pLDH) Antigen Test		Sensitivity %	Specificity %	PPV (%)	NPV (%)
		Positive	Negative				
Malaria Pf Positive Whole Blood Samples	19	19	0	100	-	100	-
Malaria Pv Positive Whole Blood Samples	31	31	0	100	-	100	-
Malaria Negative Whole Blood Samples	150	0	150	-	100	-	100

LIMITATIONS :

- As with all diagnostic tests, the test result must always be correlated with clinical findings.
- The results of test are to be interpreted within the epidemiological, clinical and therapeutic context. When it seems indicated, the parasitological techniques of reference should be considered (microscopic examination of the thick smear and thin blood films).
- Any modification to the above procedure and / or use of other reagents will invalidate the test procedure.
- Interference due to presence of heterophile antibodies in patient's sample can lead to erroneous analyte detection in immunoassay, has been reported in various studies. uses HETEROFILIC BLOCKING REAGENT (HBR) to inhibit majority of these interferences.
- In case of mixed infection (*P. falciparum* with other malarial species), both, 'Pf' and 'Pan' malaria band will be positive. Hence, differentiation of infection due to *P. vivax*, *P. ovale* or *P. malariae* cannot be done.
- While monitoring therapy, using the 'Pan' band, if the reaction of the test remains positive with the same intensity after 5-10 days, post treatment, the possibility of a resistant strain of malaria has to be considered.
- In *P. falciparum* malaria infection, HRP-2 is not secreted in the gametogamy stage. Hence, in "Carriers", the HRP-2 band may be absent.
- HRP-2 levels, post treatment persists up to 15 days, the 'Pan' band can be used to monitor success of therapy, in *P. falciparum* malaria cases.
- In a few cases, where the HRP-2 band is positive and the 'Pan' malaria band is negative, it may indicate a case of post treatment malaria. However, such a reaction pattern may also be obtained in a few cases of untreated malaria. Retesting after 2 days is advised, in such cases.
- Do not interpret the test results beyond 30 minutes.

REFERENCES

- Howard, R.J., et al., 1986: Secretion of a Malarial Histidine-rich Protein (Pf. HRP-II) from *Plasmodium falciparum*-infected Erythrocytes. J. Cell Biol., 103, 1269-1277
- Rock, E.P., et al., 1987: Comparative Analysis of the *Plasmodium falciparum* Histidine-Rich Proteins HRP-I, HRP-II and HRP-III in Malaria Parasites of Diverse Origin. Parasitol., 95, 209-227.
- Hunte-Cooke A., et al., (1999) Comparison of a Parasite Lactate Dehydrogenase-based Immunochromatographic © Antigen Detection assay (OptiMAL) with Microscopy for the Detection of Malaria Parasites in Human Blood Samples. Am J. Trop Med 60(2). 173-176

IVD	In Vitro Diagnostic Use
	Manufacturer
	Manufacturing Date
	Expiry Date
LOT	Lot Number
	Store at 4°C to 40°C
	Single Use
	Number of tests in the pack
	Do not use if pouch or kit damaged
	This side Up
	Read package insert before use



MANUFACTURED BY

ImmunoScience India Private Limited
Gat No. 41, Kusgaon, Shivapur-Velhe Road,
Tal- Bhor, Pune, Maharashtra (India) -412205.